

a measurable difference...



Cutting Press, Pneumatic

Model: C0050

The IDM Cutting Press Pneumatic is fast and easy to operate, used for sample preparation of various different sized test pieces. Tool Steel and Ruler Dies of standard and custom sizes can be used and interchanged with ease, allowing a large range of dies to be used with the Cutting Press. Cutting dies can be fixed onto the moving top plate, or can sit directly onto the sample, allowing the user to change the dies with ease. The IDM Cutting Press Pneumatic is one of the easiest and versatile cutting presses to operate on today's market. It is compact in size and fits onto a work table to conserve work space. The Cutting Press is great for laboratory use and low volume production.

Applications:

- Paper
- Cartonboard
- Polyurethane
- Plastics
- Foil
- Textiles
- Rubber
- and many other materials

Features:

- Cutting area: 300mm x 300mm
- 60mm Opening
- 95mm clearance from press plate to cutting board
- Front load for easy operation
- One or more sheets can be cut simultaneously
- Replaceable Cutting Boards

Options:

- Custom cutter sizes available upon request
- Safety guard
- Various sized Cutting Dies can be used with the Cutting Press

Connections:

- Air supply: 80 — 100 psi

Dimensions:


- **H:** 705mm
- **W:** 370mm
- **D:** 450mm
- **Weight:** 140kg



Can be used with the following Cutting Dies



C0024

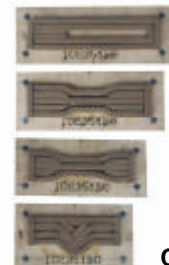

C0025

10-11 Colrado Court, Hallam, Victoria 3803 Australia.

Phone: +61 3 9708 6885

Fax: +61 3 9708 6770

Website: www.idminstruments.com.au

Email: idm@idminstruments.com.au

All sizes and dimensions shown on this data sheet are averages only and images are as accurate as the medium allows. Products and details shown may change without notice as IDM Instruments Pty Ltd[®] reserves the right to make alterations to its product range from time to time.

Copyright © 2006 IDM Instruments Pty Ltd[®]. All Rights Reserved.