

a measurable difference...



Extensiometer For Fabric

Model: E0004

The methods of tests in this standard are applicable to woven, and to warp-knitted fabrics. They are particularly applicable to stretch fabrics such as are obtained by the use of elastomeric fibres or bulked yarns, or by a process such as slack mercerisation.

A specimen of standard dimensions is stretched by applying a specified load. The extension is observed and the tension removed. The length of the unstretched specimen is observed after recovery.

Applications:

- Woven Fabric
- Knitted Fabric

Features:

- Cutting templates
- Adjustable clamps

Standards:

- BS 4952

Options:


- Spare weights

Connections:

- Manual Operation

Dimensions:

- **H:** 100mm
- **W:** 125mm
- **D:** 435mm
- **Weight:** 15kg

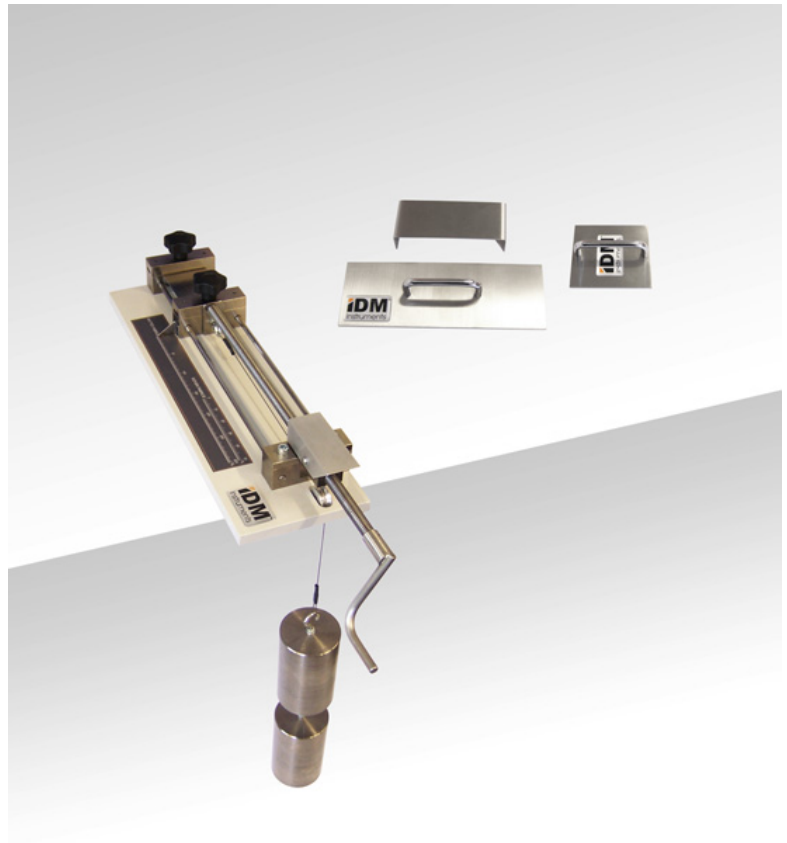

10-11 Colrado Court, Hallam, Victoria 3803 Australia.

Phone: +61 3 9708 6885

Fax: +61 3 9708 6770

Website: www.idminstruments.com.au

Email: idm@idminstruments.com.au

All sizes and dimensions shown on this data sheet are averages only and images are as accurate as the medium allows. Products and details shown may change without notice as IDM Instruments Pty Ltd[®] reserves the right to make alterations to its product range from time to time.

Copyright © 2006 IDM Instruments Pty Ltd[®]. All Rights Reserved.