

a measurable difference...



Fiber Classifier Brecht & Holl Type

Model: BH-6/12

The fiber classifier is used for the separation of pulps into fines, shives and useful fibers fractions.

The size of pulp fibers influences the mechanical properties and in the quality of paper and board.

Its design is based on a pulsating membrane, sieves and continuous water flow, and therefore the separation of the different fractions is done smoothly without any damage to the fibers or change in its original characteristics.



MAJOR TECHNICAL CHARACTERISTICS


Ring mounted stainless plate with (0,2x20) mm slots for classification of shives.

Ring mounted phosphor bronze wire screen gauze with aperture #40/in for classification of long fibers .

Ring mounted phosphor bronze wire screen gauze with aperture #120/in for classification of short fibers .

Driven by built-in electric motor and eccentric with adjustable amplitude of the pulsating membrane for 6 or 12 mm

Parts in contact with water and pulp made of corrosion resistant materials.

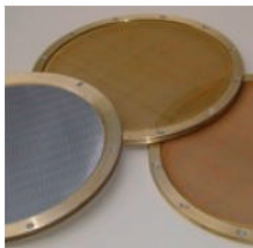

Optional accessories:
- Wire tensioning kit for easy and accurate exchange of worn screens
Applicable Standard:
Merkblatt Nr. VI/1/66

Power supply	220 VCA, 50/60 Hz, 900 W
Water supply	0,5 bar, clean, pure
Dimensions	(600 x 500 x 500) mm (W x L x H)
Weight	50 kg

Note: due to constant development our equipment design and specifications are subject to change without notice.

10-11 Colrado Court, Hallam, Victoria 3803 Australia.

Phone: +61 3 9708 6885

Fax: +61 3 9708 6770

Website: www.idminstruments.com.au

Email: idm@idminstruments.com.au

All sizes and dimensions shown on this data sheet are averages only and images are as accurate as the medium allows. Products and details shown may change without notice as IDM Instruments Pty Ltd[®] reserves the right to make alterations to its product range from time to time.

Copyright © 2006 IDM Instruments Pty Ltd[®]. All Rights Reserved.